

RESEARCH

Open Access

The application of multiple reaction monitoring and multi-analyte profiling to HDL proteins

Hussein N Yassine^{1*}, Angela M Jackson², Chad R Borges³, Dean Billheimer⁴, Hyunwook Koh¹, Derek Smith³, Peter Reaven⁵, Serrine S Lau⁶ and Christoph H Borchers^{2,7}

Abstract

Background: HDL carries a rich protein cargo and examining HDL protein composition promises to improve our understanding of its functions. Conventional mass spectrometry methods can be lengthy and difficult to extend to large populations. In addition, without prior enrichment of the sample, the ability of these methods to detect low abundance proteins is limited. Our objective was to develop a high-throughput approach to examine HDL protein composition applicable to diabetes and cardiovascular disease (CVD).

Methods: We optimized two multiplexed assays to examine HDL proteins using a quantitative immunoassay (Multi-Analyte Profiling- MAP) and mass spectrometric-based quantitative proteomics (Multiple Reaction Monitoring-MRM). We screened HDL proteins using human xMAP (90 protein panel) and MRM (56 protein panel). We extended the application of these two methods to HDL isolated from a group of participants with diabetes and prior cardiovascular events and a group of non-diabetic controls.

Results: We were able to quantitate 69 HDL proteins using MAP and 32 proteins using MRM. For several common proteins, the use of MRM and MAP was highly correlated ($p < 0.01$). Using MAP, several low abundance proteins implicated in atherosclerosis and inflammation were found on HDL. On the other hand, MRM allowed the examination of several HDL proteins not available by MAP.

Conclusions: MAP and MRM offer a sensitive and high-throughput approach to examine changes in HDL proteins in diabetes and CVD. This approach can be used to measure the presented HDL proteins in large clinical studies.

Keywords: High density lipoprotein, Proteomics, Multiple reaction monitoring, Multi-analyte panel, Diabetes, Cardiovascular disease

Introduction

Recent findings suggest that HDL carries over 80 proteins involved in lipid metabolism, complement regulation, acute phase response and proteinase inhibition [1]. This protein cargo appears to be remodeled in metabolic syndrome [2], cardiovascular disease [1,3] and after cholesterol therapies [4]. However, measurement of proteins associated with lipids and that span several orders of magnitude in abundance can be challenging. Conventional LC-MS/MS assays and label-free quantitation methods using spectral counting or extracted ion chromatograms (XIC) used in the previous studies [1-4] can be lengthy

and are limited to small sample sizes. Multiple Reaction Monitoring (MRM), on the other hand, is a tandem MS (MS/MS) scan mode unique to triple quadrupole MS instrumentation that is capable of rapid, sensitive, and specific quantitation of peptides in highly complex sample matrices, such as plasma [5,6]. MRM is a targeted approach that requires knowledge of the molecular weight the peptide of interest and its fragmentation pattern, leading to the generation of target “transitions” for monitoring protein levels. When coupled with stable isotope peptide standards (SIS peptides), quantitation using MRM can be highly reproducible [7]. MRM quantitation has been successfully applied to plasma proteins [5,6] and, more recently, to a limited subset of HDL proteins [8]. Without pre analytical depletion or enrichment, MRM sensitivity can be compromised by the

* Correspondence: hyassine@usc.edu

¹Department of Medicine, University of Southern California, Los Angeles, CA, USA

Full list of author information is available at the end of the article

existence of high abundance proteins. Although the sensitivity of MRM allows analyses at levels down to one ng/mL, the accuracy at the lower end of the abundance range can be problematic. In this case, the use of methods that deplete the higher abundance proteins is often needed. These methods, however, can lead to the unintentional removal of proteins that are attached to the depleted proteins. An alternative approach is the use of immunoassays that offer much improved sensitivity. One example of these assays is multiplexed immunoassay panels (MAP). MAP complements MRM at this low end and helps ensure that proteins are fully and correctly quantitated at the low range.

There has been considerable interest in understanding HDL functions in light of the strong inverse association of HDL cholesterol and cardiovascular disease (CVD) and the lack of improved CVD outcomes after raising HDL's cholesterol content in three recent trials [9-11]. Knowledge of HDL protein composition promises to improve our understanding of HDL's multiple functions in CVD. In this study, we optimized two approaches to measure HDL proteins using multiplexed immunoassays (Multi-Analyte Profiling-MAP) and quantitative proteomics (MRM) platforms. We extended the application of these two techniques to HDL isolated from participants with and without diabetes and CVD. Our findings suggest that MAP can be used to monitor low abundance proteins on HDL, whereas MRM allows examining higher abundance HDL proteins without the need for pre-existing antibodies.

Materials and methods

Clinical samples

The study was approved by the University of Arizona Institutional Review Board, and all patients provided written informed consent prior to testing. The samples were collected at University of Arizona Medicine Clinics. Participants reported to the Center for Clinical and Translational Sciences (CaTS) after an overnight fast. Samples were collected in EDTA containing tubes. The presence of cardiovascular disease was defined by a prior history of coronary artery bypass surgery (CABG), percutaneous transluminal angioplasty (PTCA), prior MI, or thrombotic stroke. Patients were older than 18 years of age. New diagnosis of diabetes was based on 2 hour oral glucose tolerance test (OGTT) > 200 mg/dl or glycated hemoglobin (HbA1c > 6.5%). Established diabetes was defined by clinical history. The study included disease free participants with advanced kidney disease, diabetes and CVD. The study excluded subjects if they met any of the following criteria: had type 1 diabetes, were on an active weight loss program, history of cancer, HIV, or steroid use.

HDL Isolation technique

HDL isolation by centrifugation was based on a modification of a previously published protocol [1]. In brief, potassium bromide (KBr ~55 mg) was added to 310 μ L of plasma samples to create a density of 1.21 g/mL. The sample was overlaid with 200 μ L of 1.21 g/mL density solution for a total volume was 500 μ L. Samples were then spun at 120,000 rpm, at 16°C for 2 hours (Beckman TLX ultracentrifuge with a type 120.1 fixed angle rotor using thick-walled 500 μ L Polycarbonate tubes, item 343776). The upper 125 μ L solution that had a density of less than 1.21 g/mL was removed and 150 μ L of NaCl/EDTA solution (0.9% (w/v) NaCl, 0.1% (w/v) EDTA, pH 7.4) was added to each tube for a final density of 1.063 g/mL. Subsequently, 225 μ L of 1.06 KBr solutions in NaCl/EDTA was underlaid with a final volume of 500 μ L for a second 2 hour spin at the same parameters listed. The bottom 125 μ L (HDL fraction) of solution was removed for further analysis. Four HDL samples were sent to Myriad RBM to externally validate our measurements using an immunoassay in a CLIA certified laboratory. To confirm depletion of albumin and apoB proteins from HDL samples, HDL from a CVD pool was isolated with a second approach that involved long centrifugations at two sequential spins each of 10 hours duration using the above technique. Albumin levels were then measured using a commercial ELISA (Assaypro).

Measurements of HDL proteins

We screened HDL using the human MAP panel (90 proteins) and MRM panel (56 proteins). The proteins that were detected on HDL are summarized in Table 1. Four HDL samples were used to compare protein measures by MAP and MRM. Subsequently, two pools of HDL samples (each pool 500 μ L combined from 10 HDL isolates) - one pooled from ten non-diabetic subjects defined as the control group and a second pooled from ten subjects with both diabetes and CVD defined as the disease group were run on the MAP platform using the HumanMAP panel and by MRM.

MAP

The samples were analyzed at Myriad RBM that uses Luminex xMAP. Luminex xMAP is a well-established particle array system that is based on beads with unique fluorescent signatures with proprietary multi-analyte panel targets assessed for cross-reactivity. This technology has been used for the detection of numerous targets, such as cytokines [12], cancer markers [13], and indicators for various disease states [14]. The assay details of this platform are well described in Myriad RBM website (www.myriadrbm.com). Here we applied MAP to HDL samples obtained by ultracentrifugation in both 2X2 and 10X10 isolation techniques, calculating the limit of detection and the recovery analysis after 12 dilutions. The data is

Table 1 Proteins that were quantified on HDL from the healthy and diseased sample pool

MAP (69 proteins quantified out of 90)		MRM (32 proteins quantified out of 56)	
Protein	Uniprot ID	ID	Uniprot ID
C-Reactive Protein (CRP)	P02741	Apolipoprotein D	P05090
von Willebrand Factor (vWF)	P04275	Apolipoprotein A-I	P02647
Adiponectin	Q15848	Apolipoprotein A-II precursor	P02652
Fibrinogen	P02671, P02675, P02679	Apolipoprotein A-IV	P06727
Serum Amyloid P-Component (SAP)	P02743	Apolipoprotein B-100	P04114
Sex Hormone-Binding Globulin (SHBG)	P04278	Apolipoprotein C-I lipoprotein	P02654
Immunoglobulin A (IgA)	P44969	Apolipoprotein C-III	P02656
Complement C3 (C3)	P01024	Apolipoprotein E	P02649
Protein S	P26447, P06703	Apolipoprotein L1	O14791
Thrombospondin-1	P07996	Beta-2-glycoprotein I_Apo H	P02749
Haptoglobin	P00738	Apolipoprotein C-II	P02655
Complement factor H	P08603	Apolipoprotein(a)	P08519
Ferritin (FRTN)	P02794, P02792	Apolipoprotein M	O95445
Alpha-2-Macroglobulin (A2Macro)	P01023	PON 1	P27169
Myeloperoxidase (MPO)	P05164	Albumin	P02768
Plasminogen Activator Inhibitor 1 (PAI-1)	P05121	Alpha-1-Anti-trypsin	P01009
Tumor Necrosis Factor Receptor 2 (TNFR2)	P20333	Alpha-1B-glycoprotein	P04217
Insulin	P01308	Alpha-2-antiplasmin	P08697
Vitronectin	P04004	Alpha-2-HS-glycoprotein	P02765
Calcitonin	P01258	Clusterin	P10909
Beta-2-Microglobulin (B2M)	P61769	Complement C1 inactivator	P05155
CD5L	O43866	Complement C3	P01024
Intercellular Adhesion Molecule 1 (ICAM-1)	P05362	Complement C4 beta chain	P0C0L5
CD 40 antigen (CD40)	Q6P2H9	Complement C4 gamma chain	P0C0L5
Carcinoembryonic Antigen (CEA)	P06731	Complement C9	P02748
Vascular Endothelial Growth Factor (VEGF)	P15692	Complement factor H	P08603
Thyroxine-Binding Globulin (TBG)	P05543	Fibrinogen alpha chain	P02671
Vascular Cell Adhesion Molecule-1 (VCAM-1)	P19320	Fibrinogen beta chain	P02675
Thyroid-Stimulating Hormone (TSH)	P01215, P01222	Fibrinogen gamma chain	P02679
Brain-Derived Neurotrophic Factor (BDNF)	P23560	Fibrinopeptide A	P02671
Matrix Metalloproteinase-3 (MMP-3)	P08254	Haptoglobin beta chain	P00738
Fatty Acid-Binding Protein, heart (FABP, heart)	P05413	Hemopexin	P02790
Tissue Inhibitor of Metalloproteinases 1 (TIMP-1)	P01033	Heparin cofactor II	P05546
Myoglobin	P02144	Kininogen-1	P01042
Immunoglobulin M (IGM)	P01871	L-selectin	P14151
Interleukin-8 (IL-8)	P10145		
Interleukin-1 beta (IL-1 beta)	P01584		
EN-RAGE	P80511		
Interleukin-1 alpha (IL-1 alpha)	P01583		
Monocyte Chemotactic Protein 1 (MCP-1)	P13500		
Macrophage Inflammatory Protein-1 beta (MIP-1 beta)	P13236		
Alpha-1-Antitrypsin (AAT)	P01009		
Fetuin A (alpha HS-glycoprotein)	P02765		

Table 1 Proteins that were quantified on HDL from the healthy and diseased sample pool (Continued)

Interleukin-15 (IL-15)	P40933
Epidermal Growth Factor (EGF)	P01133
Apolipoprotein D (Apo D)	P05090
Eotaxin-1	P51671
Apolipoprotein(a) (Lp(a))	P08519
Macrophage-Derived Chemokine (MDC)	O00626
Clusterin (CLU)	P10909
Prostate-Specific Antigen, Free (PSA-f)	P07288
Apolipoprotein A-I (Apo A-I)	P02647
Leptin	P41159
Matrix Metalloproteinase-9 (MMP-9)	P14780
Interleukin-10 (IL-10)	P22301
Interleukin-18 (IL-18)	Q14116
Interleukin-2 (IL-2)	P60568
Lymphotactin	P47992
T-Cell-Specific Protein RANTES (RANTES)	P13501
Erythropoietin (EPO)	P01588
Serum Glutamic Oxaloacetic Transaminase (SGOT)	P17174
Apolipoprotein E (Apo E)	P02649
Epithelial-Derived Neutrophil-Activating Protein 78 (ENA-78)	P42830
Stem Cell Factor (SCF)	P21583
Interleukin-13 (IL-13)	P35225
Apolipoprotein B (Apo B)	P04114
Growth Hormone (GH)	P01241
Interleukin-4 (IL-4)	P05112
Interleukin-1 receptor antagonist (IL-1ra)	P18510

presented in Additional file 1: Table S1. MAP was then applied to the control and disease HDL sample pool.

MRM

The samples were analyzed at the University of Victoria - Genome BC Proteomics Centre with a dedicated core service for MRMs with the capacity of high throughput stable isotope peptide production for absolute quantification. We screened HDL protein using existing published transitions that were previously applied to plasma proteins [6] or to newly developed transitions from proteins that are associated with HDL (LCAT, CETP, PLTP, PON1, Apolipoprotein D, Apolipoprotein M, Apolipoprotein L1, Apolipoprotein CII) or of interest to CVD (Macrophage migration inhibitory factor, Matrix Gla protein). The transition list used is provided in Additional file 1: Table S2. For the generation of CV data, the samples were injected four times per method, with two methods in total. The total number of transitions per sample was 142 (method 1-low

abundant specific targets) + 88 (method 2-general high abundance plasma protein targets). There were between 1–5 transitions monitored for each peptide. Proteins with new transitions, or those that were low in abundance were monitored using either multiple peptides (as CETP, MIF, PON1) or up to 5 transitions per peptide. The methods were not scheduled as the retention times shifted between HDL isolation procedures that can result in non-quantifiable data. Based on these replicate runs, we selected one transition for quantitation analysis based on the lowest coefficient of variation by the 4 replicate “technical” runs. These transitions are summarized in Table 2. The selection process of these transitions (to exclude interferences or amino acid modifications such as phosphorylation and glycosylation) was previously described [6,15] in detail and included in the (Additional file 1: MRM methods). Representative chromatograms of the healthy pool HDL transitions and disease pool HDL transitions are also included in the figure Additional file 1: Figure S1.

Table 2 HDL peptides selected for MRM

Mass info (Q1/Q3)	Protein	Peptide sequence	Fragment ion	%CV (n = 8)
575.3/937.5	Albumin	LVNEVTEFAK	y8	5.3
555.8/797.4	Alpha-1-Anti-trypsin	LSITGTYDLK	y7	23.9
656.8/771.4	Alpha-2-antiplasmin	LGNQEPGGQTALK	y8	116.2
399.5/490.3	Alpha-2-HS-glycoprotein	HTLNQIDEVK	y4	33.6
437.2/540.3	Antithrombin-III	DDLIVSDAFHK	y++9	63.2
409.3/599.4	Apolipoprotein M	AFLTPR	y5	6.4
405.9/572.8	Apolipoprotein A-I	ATEHLSTLSEK	y10++	1.2
486.8/443.2	Apolipoprotein A-II precursor	SPELQAEAK	y++8	4.5
524.3/450.8	Apolipoprotein B-100	FPEVDLTK	y++8	20.7
516.8/466.2	Apolipoprotein C-I lipoprotein	TPDVSSALDK	y++9	5.4
519.3/666.3	Apolipoprotein C-II	TAAQNLYEK	y5	3.4
598.8/244.1	Apolipoprotein C-III	GWVTDGFSSLK	b2	14.2
436.3/659.3	Apolipoprotein D	VLNQELR	y5	6.5
484.8/588.3	Apolipoprotein E	LGPLVEQGR	y5	2.1
815.9/651.3	Apolipoprotein L1	VTEPISAESGEQVER	y++12	26.7
786.5/535.3	Apolipoprotein (a)	LFLEPTQADIALLK	y++10	45.1
511.8/751.4	Beta-2-glycoprotein I_Apo H	ATWYQGER	y6	42.6
644.8/602.3	Clusterin	ELDESQVAER	y5	45.6
501.8/731.4	Complement C3	TGLQEVEVK	y6	50.8
557.8/629.4	Complement C4 beta chain	VDGTLNLNLR	y5	89.9
362.9/487.3	Complement C4 gamma chain	ITQVLHFTK	y++8	16.1
508.6/494.3	Complement C9	TEHYEEQIEAFK	y4	44.3
671.4/830.4	Complement factor H	SPDVINGSPISQK	y8	63.0
570.8/867.5	Fibrinogen alpha chain	GSESGIFTNTK	y8	50.1
497.9/600.3	Fibrinogen gamma chain	YEASILTHDSSIR	y++11	17.3
768.8/1077.5	Fibrinopeptide A	ADSGEGDFLAEGGGVR	y11	95.1
490.8/562.3	Haptoglobin beta chain	VGYVSGWGR	y5	59.9
610.8/480.3	Hemopexin	NFPSPVDAEFR	y++9	147.9
514.8/814.4	Heparin cofactor II	TLEAQLTPR	y7	31.9
626.3/1051.5	Kininogen-1	TVGSDTFYFSK	y9	28.9
497.8/794.4	L-selectin	AEIEYLEK	y6	64.4
592.8/943.5	PON 1	IQNILTEEPK	y8	51.0

HDL proteins that were quantified using MRM. The CV was calculated using 4 technical runs.

Statistical analysis

We used the statistical program R2.1 (R core development team). The characteristics of the two pools were compared using an independent sample test or a proportion test. The correlation coefficient and the p value generated between plasma proteins and HDL proteins was obtained using spearman correlation test.

Results and discussion

Description of study participants

The goal of this project was to use sensitive and high throughout approach to analyze HDL proteins in health

and vascular disease such as diabetes and CVD. The samples were selected to detect differences in HDL proteins between healthy controls and patients with chronic vascular disease. The study demographics and biochemical measurements are listed in Table 3. The two groups were age and sex matched. All diseased individuals had a history of CVD events prior to participation. As shown in Table 3, diseased subjects were obese, with lower HDL cholesterol, elevated triglyceride levels, uncontrolled diabetes, elevated inflammation (CRP), and evidence of chronic kidney disease as revealed by the elevated plasma creatinine.

Table 3 Demographic and biochemical characteristics of study participants

	Control (n = 10)	Disease (n = 10)	p value
Age (yrs)	58.4 (6.1)	62.2 (7.5)	0.233
Sex (M:F)	4:6	5:5	0.99
BMI (kg/ m ²)	24.3 (4.5)	37.8 (7.7)	0.001
Systolic BP (mm Hg)	122.4 (9.4)	131.4 (17.4)	0.173
Diastolic BP (mmHg)	74.8 (5.7)	70.7 (9.8)	0.275
LDL (mg/dL)	138.5 (35)	114.1 (42.2)	0.178
HDL (mg/dL)	58.7 (9.5)	36.1 (7.46)	<0.001
CRP (mg/dL)	2.6 (2.27)	13.3 (11)	0.014
Triglycerides (mg/dL)	111.1 (55)	309.2 (205)	0.014
Creatinine (mg/dL)	0.7 (0.15)	1.65 (0.93)	0.011
Glycated Hemoglobin %	5.4 (0.38)	9.2 (2.71)	0.002
Urine Microalbumin (mcg/mg creatinine)	8.4 (6.7)	1667.5 (2596)	0.114

The samples were pooled from the control and disease groups. Values are means (SD).


HDL and Plasma proteins

One challenge in HDL proteomics is differentiating whether proteins are actually present within HDL or are predominately carried within plasma and have become loosely associated with HDL before or during the isolation procedure. We initially isolated HDL using a two sequential spins (2 hrs × 2) and we are able to detect several low abundance proteins commonly associated with atherogenesis and inflammation using MAP (Table 1) of participants with both diabetes and CVD. To confirm whether these proteins are part of the HDL fraction or loosely attached plasma proteins, a longer centrifugation process was employed to isolate HDL (10 hrs × 2). Plasma albumin levels were determined to be 45 mg/mL. After the (2 hrs × 2) centrifugation, albumin levels within the HDL fraction were 0.1 mg/mL. After prolonged sequential centrifugation step (10 hrs × 2), albumin levels were below the detection limits of this assay. Proteins detectable on HDL (by both centrifugation techniques) assessed by MAP, along with plasma protein concentrations are summarized in Table 4. We then correlated the concentrations of these proteins between the different fractions. We found no correlation between plasma and HDL proteins identified in either the (2 hrs × 2) or (10 hrs × 2) HDL isolation preparations. For example, Apo A-I concentrations of 0.36, 0.27, 0.18 mg/mL in plasma HDL (2 hrs × 2) and HDL (10 hrs × 2), respectively shows clear retention on HDL compared to fibrinogen with concentrations of 4.7, 0.0021, 0.0012 mg/mL. These findings confirm that non-HDL proteins were efficiently depleted with the longer centrifugation. In contrast, the HDL proteins isolated with the (2 hrs × 2) and (10 hrs × 10) methods were highly correlated (Figure 1, $r = 0.95$, $p < 0.001$). As

Table 4 Concentrations of Plasma and HDL proteins from the pooled sample of the diseased subject (mg/mL)

Protein	Plasma	HDL 2X2	HDL 10X10
Apolipoprotein A-I (Apo A-I)	0.36	0.259	0.178
Alpha-1-Antitrypsin (AAT)	1.5	0.0688	0.00736
Immunoglobulin A (IgA)	2.5	0.00422	0.00054
Immunoglobulin M (IGM)	2.1	0.00275	0.00117
Complement C3 (C3)	1.5	0.00211	0.000792
Fibrinogen	4.7	0.00209	0.00118
Haptoglobin	2.3	0.00191	0.000438
Alpha-2-Macroglobulin	1.9	0.000857	0.00104
Apolipoprotein (a) (Lp (a))	0.684	0.108	0.081
Fetuin-A	0.2	0.068	0.015
Complement Factor H	0.1	0.057	0.0036
Apolipoprotein D (Apo D)	0.05	0.041	0.018
Apolipoprotein C-III (Apo C-III)	0.07	0.039	0.028
Apolipoprotein B (Apo B)	1.0	0.028	0.01
Apolipoprotein H (Apo H)	0.246	0.016	0.00091
Vitronectin	0.26	0.013	0.0018
Apolipoprotein E (Apo E)	0.04	0.0074	0.0042
Clusterin (CLU)	0.37	0.0042	0.0011
Thyroxine-Binding Globulin (TBG)	0.044	0.0025	0.00033
Beta-2-Microglobulin	0.0041	0.00094	0.00033
Vitamin K-Dependent Protein S	0.05	0.00027	0.000074
Serum Amyloid P-Component	0.017	0.00015	0.000051
C-Reactive Protein (CRP)	0.015	0.00013	0.000023
CD5 (CD5L)	0.00053	0.000106	0.0000096
Vascular Cell Adhesion Molecule-1 (VCAM-1)	0.000761	0.000039	0.000011
Myoglobin	0.00013	0.000031	0.000016
Tissue Inhibitor of Metalloproteinases 1 (TIMP-1)	0.000104	0.0000097	0.0000012
Thrombospondin-1	0.00482	0.000008	0.0000036
Plasminogen Activator Inhibitor 1 (PAI-1)	0.000072	0.000004	0.0000005
Tumor Necrosis Factor Receptor 2 (TNFR2)	0.000014	0.0000017	0.0000001
T-Cell-Specific Protein RANTES (RANTES)	0.0000097	0.0000015	0.0000004
EN-RAGE	0.000067	0.0000001	< LOW >
Myeloperoxidase (MPO)	0.00183	< LOW >	0.000019
Sex Hormone-Binding Globulin (SHBG)	0.0016195	0.0000218	0.0000006

HDL 2.2 and 10.10 refer to two sequential ultra-centrifugations 2 or 10 hours each. Samples were diluted 12 times and were run once at each dilution. The reported concentrations were in the linear range of the assay (Additional file 1: Table S1). The strength of this technique is in the ability to measure these HDL proteins across a wide concentration range.

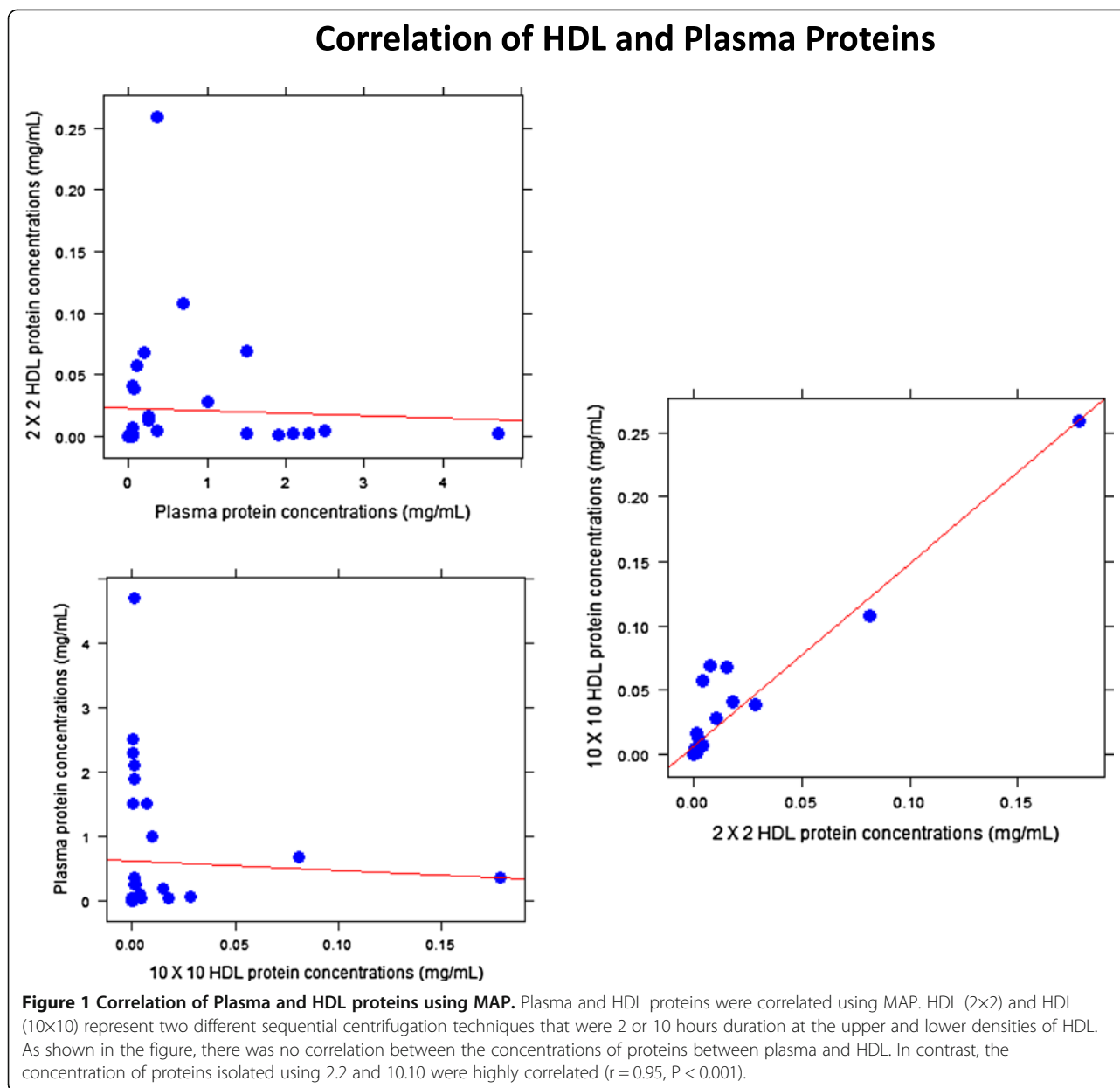

expected the concentrations of HDL proteins were greater in the (2 hrs \times 2) than the (10 hrs \times 2) technique (Table 3). The concentration of Apo A-I in plasma was low in this sample pool from diseased participants compared to healthy controls. Published Apo A -I levels are in the 0.9-1.6 mg/mL range. This may suggest significant HDL remodeling in advanced vascular disease.

Correlation of MAP and MRM

MRM's performance for low abundance proteins is limited without prior enrichment. Thus, MAP allows the measurements of low abundance proteins on HDL. There were 10 proteins common to both MAP and MRM, however, one protein (Apo E) was not detected in our samples


in 3 out of 4 samples by analyzed by MRM. We correlated 9 proteins measured by both approaches in 4 HDL samples. The results are summarized in Table 5. We did not observe a good correlation for lipoprotein (a), Apo D and Apo CIII between the two methods. Most of Lipoprotein (a) is associated with LDL and present in low abundance on HDL [16] reflecting a challenge for measuring this protein on HDL by MRM. Apo D on the other hand, has a high degree of homology to retinol binding protein or other members of the alpha 2 microglobulin protein superfamily [17]. This might potentially pose a challenge for the MAP technique. Apo CIII was another protein that was measured by MAP and MRM with a weak correlation between the two measurement platforms. It is not clear

of this approach is limited. A larger study of subjects with and without disease is needed to fully characterize the distribution of these low abundance proteins on HDL. Unlike MAP, the sensitivity of MRM is limited without prior sample enrichment. Hence, the accuracy of the assay is compromised when sample concentrations fall below the linear range of the assay as evidenced by the increased coefficient of variation (CV) in several MRM proteins. By MRM, we targeted transitions of 56 proteins. We were able to detect 32 proteins of which 11 proteins had CVs less than 20% on the replicate runs (Tables 1 and 2). Clusterin is a high abundant protein that is associated with HDL. We were able to quantitate clusterin accurately in plasma (CV < 5% [6]); however we were not able to reliably measure it in HDL as was previously demonstrated [8] reflecting the challenges of using MRMs for protein quantification when levels of a given target falls out of the dynamic range in diseased states. Comparing HDL proteins in the CVD versus control pool by MRM revealed the depletion of several lipid metabolism proteins such as Apo CI, CII and PON-1 in CVD. In contrast, the concentrations of several acute phase proteins such as clusterin, Complement C9, alpha-1-antitrypsin were increased in CVD (Figure 3).

These findings would facilitate conducting larger studies using the presented approach to examine the effect of inflammation and atherosclerosis associated proteins on HDL composition and function in individuals with diabetes and CVD.

Challenges of measuring lipid associated proteins

Measurement of lipid associated proteins has been successful using immunoassays based on nephelometry and mass spectrometry [1,4,8]. Nephelometry is not widely available and not amenable to multiplexing. Conventional mass spectrometry techniques using spectral counting or extracted ion chromatograms can be lengthy and challenging in large sample sizes. Here, we present two alternatives (MAP and MRM) that are amenable to multiplexing, and are high throughput. MAP is limited by the availability and quality of pre-existing antibody panels. MRM, on the other hand, can lose sensitivity at the lower end of abundance. Although the performance of MRM in lower abundance HDL proteins was inferior to immune based assays, the performance of MRM is likely to improve with the development of more sensitive mass spectrometry analyzers and better sample fractionation methods.

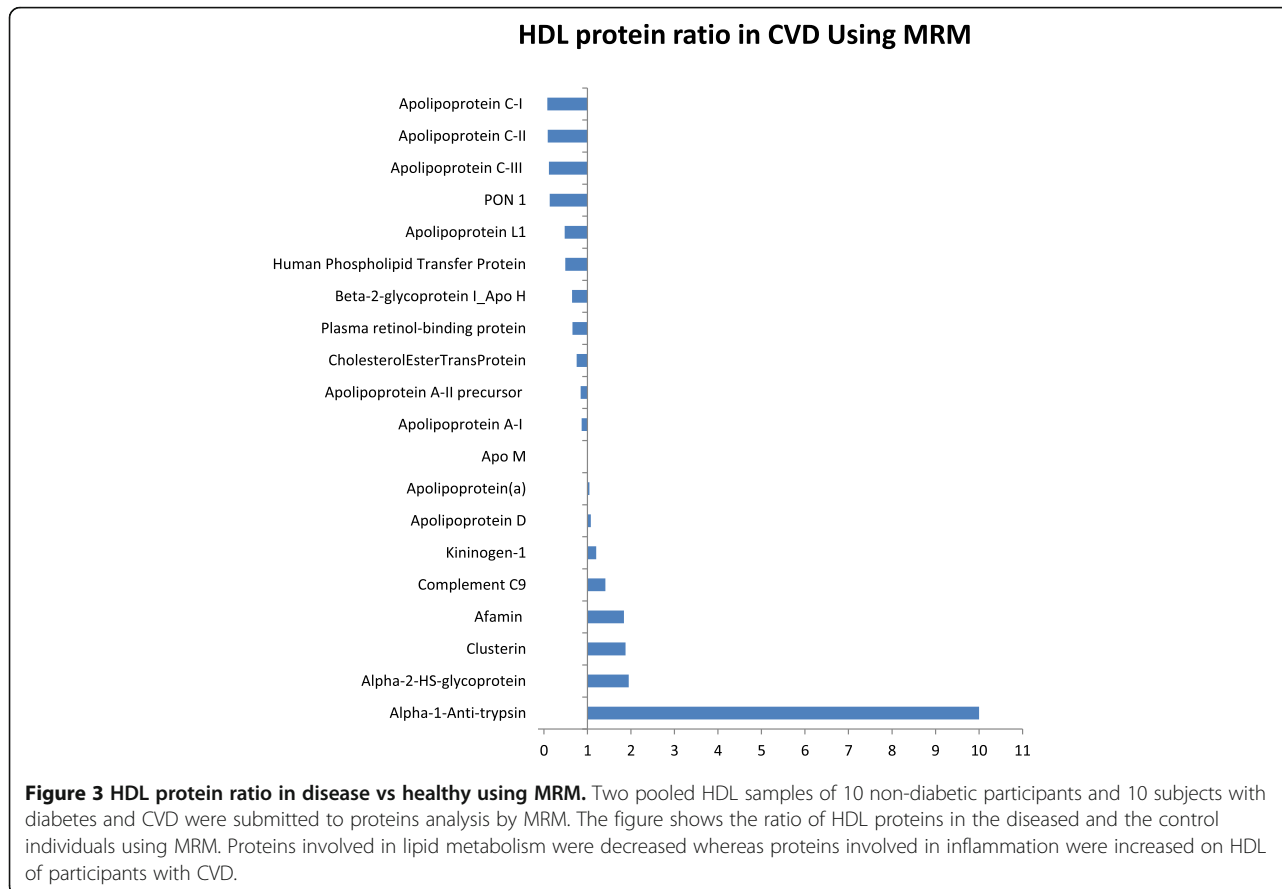

Limitations

There are several limitations of this study. First, the present study was done with two pools of clinical samples. As such, the present work describes an alternative analytical tool for the measurement of HDL proteins. Application of these assays to larger clinical data sets is feasible and will allow determination of their clinical utility. Second, the digestion efficiency of each protein monitored by MRM was not assessed. However, we have previously shown [19] that if the digestion procedure is consistent and reproducible, then the ability to compare quantitative values between samples is maintained. Thus, the results presented are better described as “relative accurate abundances”. Despite this limitation, several peptides quantitated by MRM were correlated with measurement using MAP (Table 5) affirming our claim that the MRM assay can provide absolute quantitation. A third limitation of this study is that we did not use an independent method of HDL isolation (such as reciprocal co-immunoprecipitation of a given protein with Apo A-I, or gel filtration as suggested by Davidson et al. [20]) to validate that presence of these low abundance proteins on HDL. However, HDL isolated by a longer centrifugation step had successfully depleted albumin with several of the low abundance proteins still detectable on HDL. In addition, the proteins on HDL and the corresponding plasma concentrations did not correlate, whereas the two HDL fractions were highly correlated. This finding suggests that these low abundance proteins were not contaminant plasma proteins. Our findings however, need to be replicated in a larger study group.

Conclusions

This study suggests the feasibility of measuring HDL proteins using MRM and MAP. The application of MAP and MRM to the HDL proteome offers the potential to improve our understanding of HDL functions and help direct interventions aimed at remodeling the HDL phenotype in diabetes and CVD.

Additional file

Additional file 1: lists (1) details of MRM method summary. (2) Table S1. showing the HDL protein concentrations after 12 dilutions for determination of lower limit of detection (linearity analysis) (3) Table S2. listing all the transitions used to screen for HDL proteins (4) Figure S1. showing representative chromatograms from the 4 replicate runs.

Abbreviations

HDL: High density lipoprotein; CVD: Cardiovascular disease; Apo A-I: Apolipoprotein A1; MRM: Multiple reaction monitoring; MAP: Multi-analyte profiling.

Competing interests

The author's declare that they have no competing interests.

Authors' contributions

Participated in research design: HY, SL, DS, CB, DB; conducted experiments: HY, HK, AJ; performed data analysis and interpretation: HY, DB; contributed to the writing of the manuscript: PR, CB; critically revised the manuscript: PR; all authors read and approved the final manuscript.

Acknowledgements

We would like to thank Genome Canada and Genome BC for Science and Technology Innovation Centre funding support. We thank Carol Parker for her assistance with writing the manuscript. We would also like to recognize Tyra J. Cross and Suping Zhang for the synthesis of all SIS peptides and Juncong Yang for exemplary technical support. We also thank Dr. George Tsaprailis with his assistance in running MRMs at the Arizona Proteomics Consortium. Finally, we would like to thank Jeffrey Freiser for his assistance with Myriad RBM HDL runs.

Sources of funding

Dr. Yassine was supported by K23HL107389, AHA12CRP11750017 and USC CTSI pilot UL1TR000130. The content is solely the responsibility of the authors and does not necessarily represent the official views of the National Institutes of Health. MRM mass spectrometric data were acquired by the Arizona Proteomics Consortium supported by NIEHS grant P30ES06694 to the Southwest Environmental Health Sciences Center (SWEHSC to Dr. Lau), NIH/NCI grant P30CA023074 to the Arizona Cancer Center (AZCC), and by the BIO5 Institute of the University of Arizona.

Author details

¹Department of Medicine, University of Southern California, Los Angeles, CA, USA. ²University of Victoria - Genome British Columbia Proteomics Centre, Victoria, BC, Canada. ³Arizona State University, Tempe, AZ, USA. ⁴Statistics Consulting Laboratory, University of Arizona, Tucson, AZ, USA. ⁵Phoenix VA Health Care System, Phoenix, AZ, USA. ⁶Southwest Environmental Health Sciences Center, Department of Pharmacology and Toxicology, College of Pharmacy, University of Arizona, Tucson, AZ, USA. ⁷Department of Biochemistry and Microbiology, University of Victoria, Victoria, BC, Canada.

Received: 26 September 2013 Accepted: 26 December 2013

Published: 8 January 2014

References

1. Vaisar T, Pennathur S, Green PS, Gharib SA, Hoofnagle AN, Cheung MC, Byun J, Vuletic S, Kassim S, Singh P, et al: **Shotgun proteomics implicates protease inhibition and complement activation in the antiinflammatory properties of HDL.** *J Clin Invest* 2007, **117**:746–756.
2. Hoofnagle AN, Wu M, Gosmanova AK, Becker JO, Wijsman EM, Brunzell JD, Kahn SE, Knopp RH, Lyons TJ, Heinecke JW: **Low clusterin levels in high-density lipoprotein associate with insulin resistance, obesity, and dyslipoproteinemia.** *Arterioscler Thromb Vasc Biol* 2010, **30**:2528–2534.
3. Alwaili K, Bailey D, Awan Z, Bailey SD, Ruel I, Hafiane A, Krimbou L, Laboisiere S, Genest J: **The HDL proteome in acute coronary syndromes shifts to an inflammatory profile.** *Biochim Biophys Acta* 2012, **1821**:405–415.
4. Green PS, Vaisar T, Pennathur S, Kulstad JJ, Moore AB, Marcovina S, Brunzell J, Knopp RH, Zhao XQ, Heinecke JW: **Combined statin and niacin therapy remodels the high-density lipoprotein proteome.** *Circulation* 2008, **118**:1259–1267.
5. Anderson L, Hunter CL: **Quantitative mass spectrometric multiple reaction monitoring assays for major plasma proteins.** *Mol Cell Proteomics* 2006, **5**:573–588.
6. Kuzyk MA, Smith D, Yang J, Cross TJ, Jackson AM, Hardie DB, Anderson NL, Borchers CH: **Multiple reaction monitoring-based, multiplexed, absolute quantitation of 45 proteins in human plasma.** *Mol Cell Proteomics* 2009, **8**:1860–1877.
7. Addona TA, Abbatiello SE, Schilling B, Skates SJ, Mani DR, Bunk DM, Spiegelman CH, Zimmerman LJ, Ham AJ, Keshishian H, et al: **Multi-site assessment of the precision and reproducibility of multiple reaction monitoring-based measurements of proteins in plasma.** *Nat Biotechnol* 2009, **27**:633–641.
8. Hoofnagle AN, Becker JO, Oda MN, Cavigiolio G, Mayer P, Vaisar T: **Multiple-reaction monitoring-mass spectrometric assays can accurately**

- measure the relative protein abundance in complex mixtures. *Clin Chem* 2012, **58**:777–781.
9. Boden WE, Probstfield JL, Anderson T, Chaitman BR, Desvignes-Nickens P, Koprowicz K, McBride R, Teo K, Weintraub W: **Niacin in patients with low HDL cholesterol levels receiving intensive statin therapy.** *N Engl J Med* 2011, **365**:2255–2267.
 10. Nissen SE, Tardif JC, Nicholls SJ, Revkin JH, Shear CL, Duggan WT, Ruzyllo W, Bachinsky WB, Lasala GP, Tuzcu EM: **Effect of torcetrapib on the progression of coronary atherosclerosis.** *N Engl J Med* 2007, **356**:1304–1316.
 11. Schwartz GG, Olsson AG, Abt M, Ballantyne CM, Barter PJ, Brumm J, Chaitman BR, Holme IM, Kallend D, Leiter LA, et al: **Effects of dalcetrapib in patients with a recent acute coronary syndrome.** *N Engl J Med* 2012, **367**:2089–2099.
 12. Carson RT, Vignali DA: **Simultaneous quantitation of 15 cytokines using a multiplexed flow cytometric assay.** *J Immunol Methods* 1999, **227**:41–52.
 13. Arellano-Garcia M, Hu S, Wang J, Henson B, Zhou H, Chia D, Wong D: **Multiplexed immunobead-based assay for detection of oral cancer protein biomarkers in saliva.** *Oral diseases* 2008, **14**:705–712.
 14. Morgan E, Varro R, Sepulveda H, Ember JA, Apgar J, Wilson J, Lowe L, Chen R, Shivraj L, Agadir A: **Cytometric bead array: a multiplexed assay platform with applications in various areas of biology.** *Clin Immunol* 2004, **110**:252–266.
 15. Yassine H, Borges CR, Schaab MR, Billheimer D, Stump C, Reaven P, Lau SS, Nelson R: **Mass spectrometric immunoassay and MRM as targeted MS-based quantitative approaches in biomarker development: Potential applications to cardiovascular disease and diabetes.** *Proteomics Clin Appl* 2013, **7**:528–540.
 16. Scanu AM, Fless GM: **Lipoprotein (a). Heterogeneity and biological relevance.** *J Clin Invest* 1990, **85**:1709.
 17. DRAYNA DT, McLEAN JW, WION KL, TRENT JM, DRABKIN HA, LAWN RM: **Human Apolipoprotein D Gene: Gene Sequence, Chromosome Localization, and Homology to the α 2u-Globulin Superfamily.** *Dna* 1987, **6**:199–204.
 18. Niederkofler EE, Tubbs KA, Kiernan UA, Nedelkov D, Nelson RW: **Novel mass spectrometric immunoassays for the rapid structural characterization of plasma apolipoproteins.** *J Lipid Res* 2003, **44**:630–639.
 19. Proc JL, Kuzyk MA, Hardie DB, Yang J, Smith DS, Jackson AM, Parker CE, Borchers CH: **A quantitative study of the effects of chaotropic agents, surfactants, and solvents on the digestion efficiency of human plasma proteins by trypsin.** *J Proteome Res* 2010, **9**:5422–5437.
 20. Davidson WS, Silva RA, Chantepie S, Lagor WR, Chapman MJ, Kontush A: **Proteomic analysis of defined HDL subpopulations reveals particle-specific protein clusters: relevance to antioxidative function.** *Arterioscler Thromb Vasc Biol* 2009, **29**:870–876.

doi:10.1186/1476-511X-13-8

Cite this article as: Yassine et al.: The application of multiple reaction monitoring and multi-analyte profiling to HDL proteins. *Lipids in Health and Disease* 2014 **13**:8.

Submit your next manuscript to BioMed Central and take full advantage of:

- Convenient online submission
- Thorough peer review
- No space constraints or color figure charges
- Immediate publication on acceptance
- Inclusion in PubMed, CAS, Scopus and Google Scholar
- Research which is freely available for redistribution

Submit your manuscript at
www.biomedcentral.com/submit

